



DROP SEAT

These adjustable drop seats lower the wheelchair seat so the patient's feet can touch the floor. They attach to the wheelchair seat frame using polypropylene straps and square metal rings. The seat height may be lowered from 1/4" (5/8 cm) to 3" (8 cm) depending on the wheelchair model, and/or can be adjusted to create a "tilt" or "wedge effect" to help prevent forward sliding. Made of durable, wipe-clean polyvinylchloride (PVC).

- No tools required to adjust seat height.

REF 654018 Fits most 18" (46 cm) wheelchairs

REF 654020 Fits most 20" (51 cm) wheelchairs

MONITORING

AWARNING Always follow your facility's policies and procedures for patient assessment, monitoring and rehabilitation

INSPECTION

AWARNING Inspect before each use: check for broken stitches or parts; torn, cut or frayed material; or buckles that do not hold securely. DO NOT use soiled or damaged products.

APPLICATION INSTRUCTIONS:

1. Remove the vinyl seat of the wheelchair by unscrewing the seat screws.
2. Release the polypropylene straps of the Drop Seat from the square loop rings.
3. Orient the Drop Seat to the wheelchair as shown in the photo. Wrap each end of the polypropylene straps (on each side of the chair) around the metal bar of the seat frame and thread back through the slot, pulling the ends toward the center of the Drop Seat (Fig. 1).
4. Attach the straps underneath the seat using the square loop rings (Fig. 2). All straps should pass under the Drop Seat.

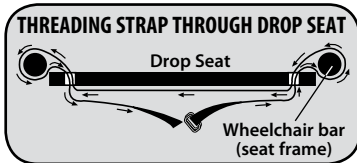


Fig. 1

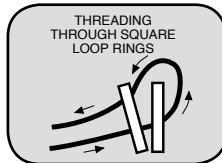


Fig. 2

CLEANING INSTRUCTIONS:

- Wipe Clean with mild detergent. OSHA approved intermediate level disinfectants can be used per manufacturer instructions. DO NOT use phenol and benzyl based disinfectants.
- After cleaning, products MUST be rinsed with water to remove any residual chemicals.
- Make sure products are completely dry before use.



WIPE CLEAN

5. Adjust the polypropylene straps so the Drop Seat is at the desired height (If desired, straps may be adjusted to create a "wedge").
6. Test the attachment for strength and security by applying weight with your foot, knee or hand down on the seat.
CAUTION DO NOT sit down until you have made sure the attachment is secure!
7. Insert a cushion on top of the Drop Seat to ensure patient comfort and proper positioning.
CAUTION DO NOT use the Drop Seat without a cushion!
NOTE: When installed and used properly, this product has a weight capacity of 250 lbs. (113 kgs.).

DROP SEAT CUSHION

This cushion is specially designed to fit the Posey 654018 Drop Seat only for 18" (46 cm) wheelchairs. Made of durable polyurethane foam and covered in blue vinyl. This cushion has notches cut out to fit over the cross bars.



REF 7020F 15½" x 16" x 2" thick (39 cm x 41 cm x 5 cm)



APPLICATION INSTRUCTIONS

1. Make sure the cushion dimensions match those of the wheelchair.
2. Place cushion on the wheelchair seating surface.
3. Align the front edge of the cushion with the front edge of the wheelchair seat.

DISPOSAL

AWARNING Properly dispose of the product per facility's policy for BIOHAZARDOUS materials.

STORAGE AND HANDLING

This device is designed for use in normal indoor environments. This device may be stored in ambient warehouse temperatures at normal humidity levels. Avoid excess moisture or high humidity that may damage product materials.