



Quality Products since 1937

PoseyBED™

POSEY HEADBOARD AND FOOTBOARD PADS


HEADBOARD AND FOOTBOARD PADS

The versatile Posey Headboard and Footboard Pads are designed to help the clinician prevent patient injury. The pads cushion the headboard and footboard, as well as helping close the gaps between the headboard and footboard and the mattress. Designed for use with the Posey Bed, they also work independently to help provide an extra measure of safety and protection for all patients. Made of durable, plush foam and covered in water-repellent, fire retardant vinyl. Pads attach to the bed using propylene straps with quick-release buckles.

REF 8002 35" H x 17" L x 2" Thick (89 cm x 43 cm x 5 cm)

APPLICATION INSTRUCTIONS:

1. Position the Headboard/footboard pads on the inside of the headboard and footboard. Ensure that the top rail of the headboard and footboard are completely covered by the extended width portion of the pad.
2. Wrap both sets of vertical straps (fig. 1) around the headboard/footboard and connect using the quick-release buckles.
3. Wrap the horizontal set of straps around the headboard/footboard and connect using the quick-release buckle.
4. Pull all straps and quick-release buckles snug, but not so tight as to disfigure the product.

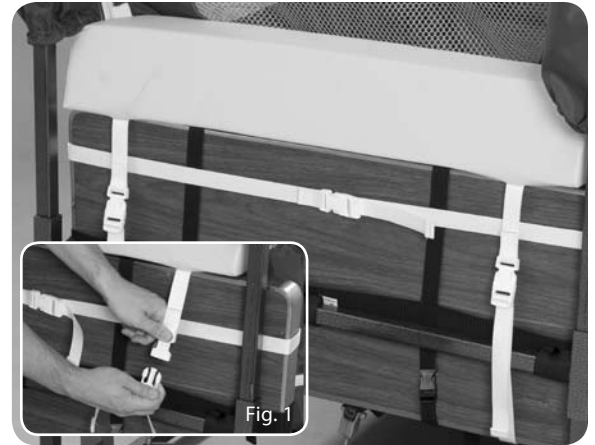

CLEANING INSTRUCTIONS:

Wipe clean with disinfectant or with warm soap and water.



Posey Company • 5635 Peck Road, Arcadia, CA 91006-0020 USA
Phone: 1.800.447.6739 • Fax: 1.800.767.3933 • Web: www.posey.com

© 2009 Posey Company. All rights reserved.

M6081 012609



Quality Products since 1937

PoseyBED™

POSEY HEADBOARD AND FOOTBOARD PADS


HEADBOARD AND FOOTBOARD PADS

The versatile Posey Headboard and Footboard Pads are designed to help the clinician prevent patient injury. The pads cushion the headboard and footboard, as well as helping close the gaps between the headboard and footboard and the mattress. Designed for use with the Posey Bed, they also work independently to help provide an extra measure of safety and protection for all patients. Made of durable, plush foam and covered in water-repellent, fire retardant vinyl. Pads attach to the bed using propylene straps with quick-release buckles.

REF 8002 35" H x 17" L x 2" Thick (89 cm x 43 cm x 5 cm)

APPLICATION INSTRUCTIONS:

1. Position the Headboard/footboard pads on the inside of the headboard and footboard. Ensure that the top rail of the headboard and footboard are completely covered by the extended width portion of the pad.
2. Wrap both sets of vertical straps (fig. 1) around the headboard/footboard and connect using the quick-release buckles.
3. Wrap the horizontal set of straps around the headboard/footboard and connect using the quick-release buckle.
4. Pull all straps and quick-release buckles snug, but not so tight as to disfigure the product.

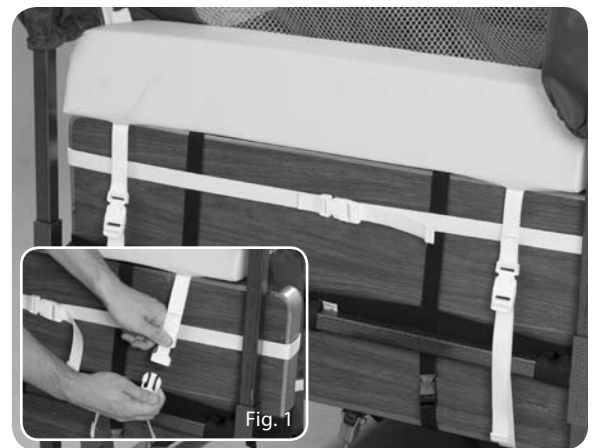

CLEANING INSTRUCTIONS:

Wipe clean with disinfectant or with warm soap and water.



Posey Company • 5635 Peck Road, Arcadia, CA 91006-0020 USA
Phone: 1.800.447.6739 • Fax: 1.800.767.3933 • Web: www.posey.com

© 2009 Posey Company. All rights reserved.

M6081 012609